




## **ProClean**

Protective hood with / without aerosol extraction

## Height adjustable protective hood

Maintenance free, pneumatically controlled, height-adjustable hood ProShield M for active user protection

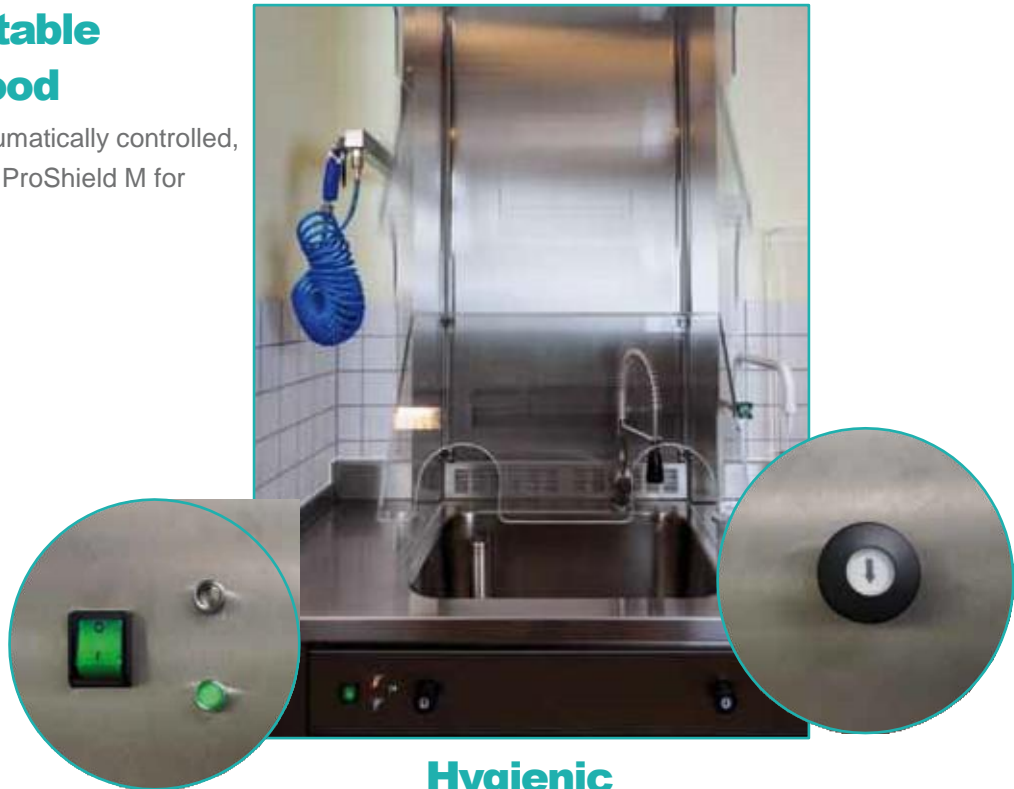

## Clear information

Clear signalization of the operating state from ProClean Air by LED

## Hygienic opening and closing

The protective hood ProShield M can be closed and opened by touch on the button. The buttons can be used either by hand or with the hip

## Measure device for disinfectant solutions

To find exact measure from disinfectant solutions



## Suction unit

ProClean Air suction unit removes the humidity from the air and removes the aerosols



## Folding protective hood

Folding protective hood ProShield S for ultrasonic washbasins



## Special hinges

Special hinges keep ProShield S safely open







### **Spray guns**

Spray guns for compressed air or water mounted on the media panel



### **Cleaning shower**

Cleaning shower with small height



### **Water tap for disinfection liquid**

Disinfection liquid resistant water tap



## ProClean

The system **Proclean** has been created for the health protection of the user during the manual cleaning procedure of instruments and endoscopes.

**Proclean** unites functionality, usefulness and hygienic comfort at the highest standard.

The concept from **Proclean** makes it possible to use the vertical moveable protective hood **Proshield M** for the cleaning basins without touching the protective hood. You can touch a footswitch or a button for opening or closing the protective hood. The germ transmission is therefore counteracted effectively


### more hygienic protection is not possible

With **Proclean air** the humid air and air with aerosols will be sucked off from the basins and will be filtered by a 2 phases filter system. The noxious materials are safely filtered out of the loaded air of the cleaning washbasin, before it is released again in the surroundings.

According the requirement, you can make the decision between **Proshield M** or **Proshield s** with or without the **Proclean air** suction unit.

The high effective combination between the protective hoods with or without the extraction units and the medical furnishings of the **Proline** series sets no limits to the equipment possibilities – from stainless steel (AISI 304) up to the electrolytically galvanized steel – everything is possible.

Numerous combination variants – from the equipment up to the accessories – ensure optimal work with maximum and effective space utilization. Practical. Uncompromising. Modern.

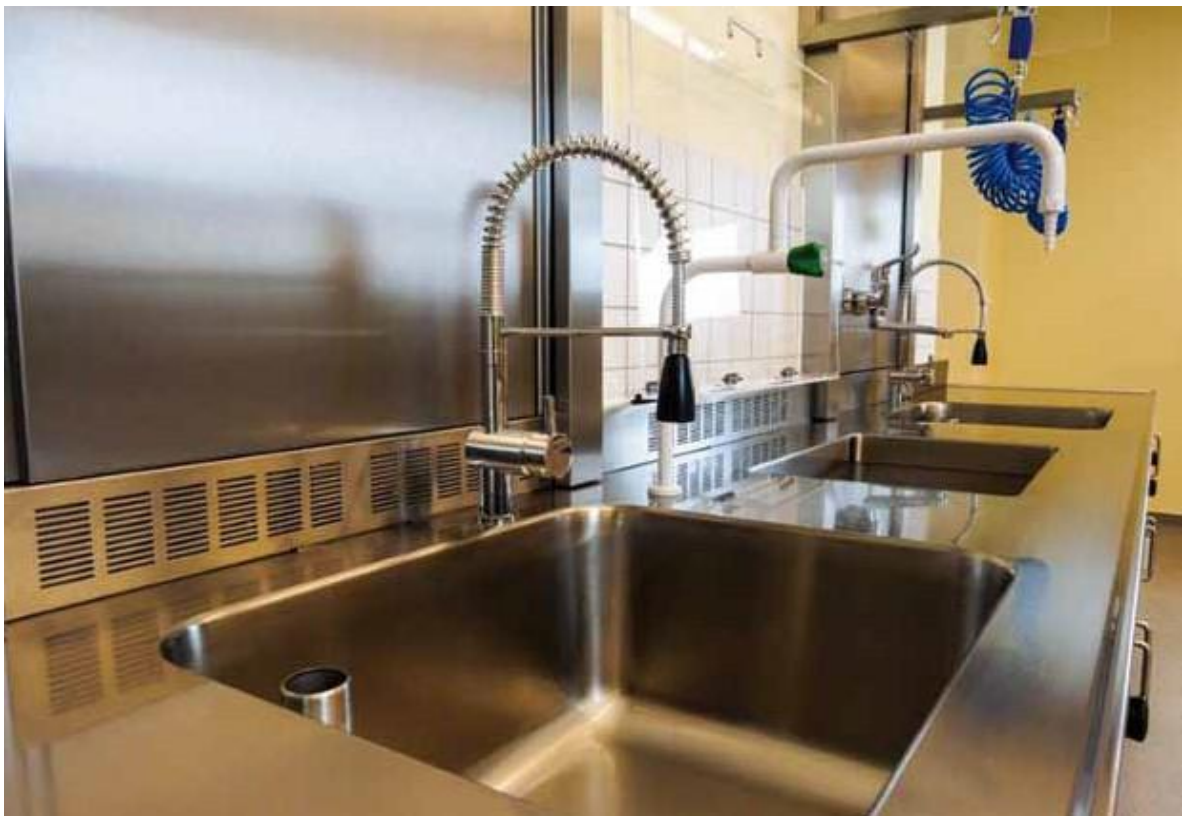

# Internally produced furnishings

- Sales
- Project management
- CAD construction
- Fabrication
- Assembly
- Service



**Lischka** GmbH

Boxberger Straße 19  
12681 Berlin  
Germany

Tel. +49 30 98 19 58-0  
Fax +49 30 98 19 58-200

info@made-by-lischka.de  
www.made-by-lischka.de



<http://made-by-lischka.de/en/schutzhauben/>

Rev.: 2/Date 08.2019

C